



IG SURGICAL

AUCKLAND PATIENT RETURN ELECTRODE

The new standard in Diathermy Safety



OKLand®



Auckland Pad

The Auckland Patient Return Electrode is a large current-limiting gel pad.

It provides increased safety to the patient and acts as a fully waterproof and reusable pressure reduction pad.

It removes the need for disposable sticky plates during diathermy procedures, removing any risk of patient harm due to pad-site burns.


The Pad is filled with a viscoelastic polymer which effectively reduces patient sores and can be used with a wide range of patient warming systems.

Auckland Pad - Dual Cord

The Dual Cord Pad has two diathermy connectors, enabling two generators to be used at the same time.

This considerably improves ease for staff, as it ends the need for multiple sticky pad placement sites.

This brings further safety benefits to the patient and additional cost savings to the hospital.

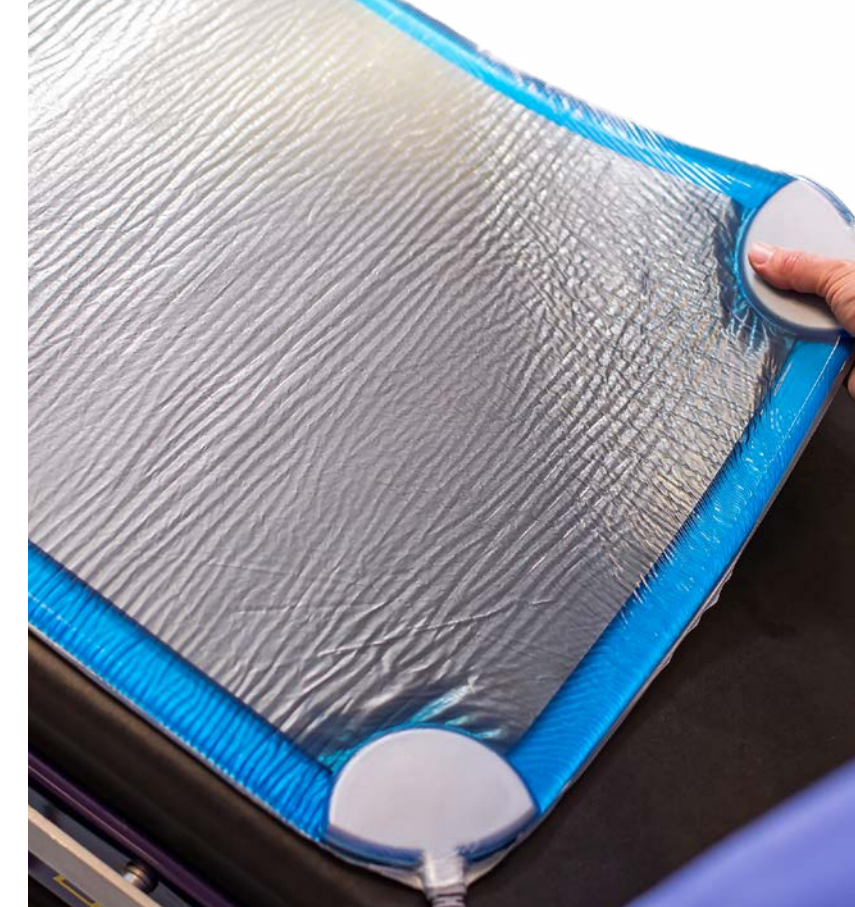

FEATURES & BENEFITS



Improved Safety: Large size of return electrode and current-limiting technology make the pad safer than sticky plates. No damage to fragile skin and it eliminates the risk of unintended burns.



Latex Free



Increased Simplicity: Easy set up with no need to find and prepare a suitable site for sticky plates. Compatible with most major brands of diathermy generators.



Radiolucent



Financial Savings: A 2 year warranty means this reusable pad can significantly reduce theatre costs including a massive reduction in clinical waste.



Full two year warranty

Auckland Pad - Paediatric

Specifically designed for babies and infants weighing from 350g. The current-limiting technology makes the safety advantages exceptional.

The ability to perform mono-polar surgery without the need to attach sticky plates onto the delicate skin of the infant, and thereby remove the risk of an unintended burn is revolutionary.



Auckland Pad - Lite

The Auckland Pad - Lite offers all the same patient safety benefits without the weight that comes with a full pressure reduction pad.

At 120cm x 50cm it is a full size pad and yet at a thickness of only 0.3cm it weighs just 1.92kg.





GET IN TOUCH!



IG SURGICAL



+44 (0)1242 262680



enquiries@igsurgical.co.uk



www.igsurgical.co.uk



IG Surgical Office:

Inter Global Surgical LLP
Suite 100
20 Winchcombe Street
Cheltenham
GL52 2LY

Description	Product Code
Patient Return Electrode Pad - Single Cord	PRE-0830
Patient Return Electrode Pad - Dual Cord	PRE-0835
Patient Return Electrode Pad - Paediatric	PRE-0840
Patient Return Electrode Pad - Lite	PRE-0845
Patient Return Electrode Pad - Lite Dual Cord	PRE-0850
Connecting Cable 4.1m (Square Head, Single Pole)	CC-01
Connecting Cable 4.1m (Square Head, Dual Pole)	CC-02
Connecting Cable 4.1m (Microphone Head, Single Pole)	CC-03
Connecting Cable 4.1m (Microphone Head, Dual Pole)	CC-04
Connecting Cable 4.1m (Flat Head, Single Pole)	CC-05
Connecting Cable 4.1m (Flat Head, Dual Pole)	CC-06
Connecting Cable 4.1m (Extended Square Head, Single Pole)	CC-07
Connecting Cable 4.1m (Extended Square Head, Dual Pole)	CC-08